
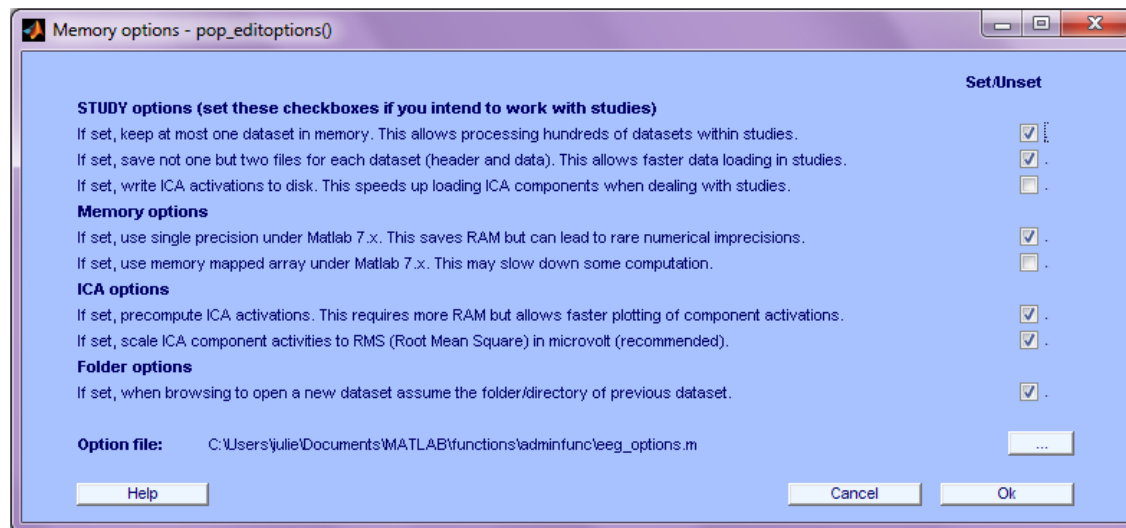
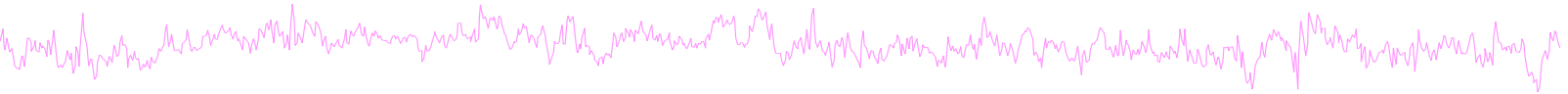


STUDY scripts

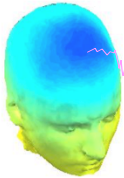




Most important option:

- Allows only one dataset to be loaded at once.
- Most STUDYs are too big to have all data loaded at once.

% Set memory options:

```
pop_edioptions( 'option_storedisk', 1);
```

Edit dataset info

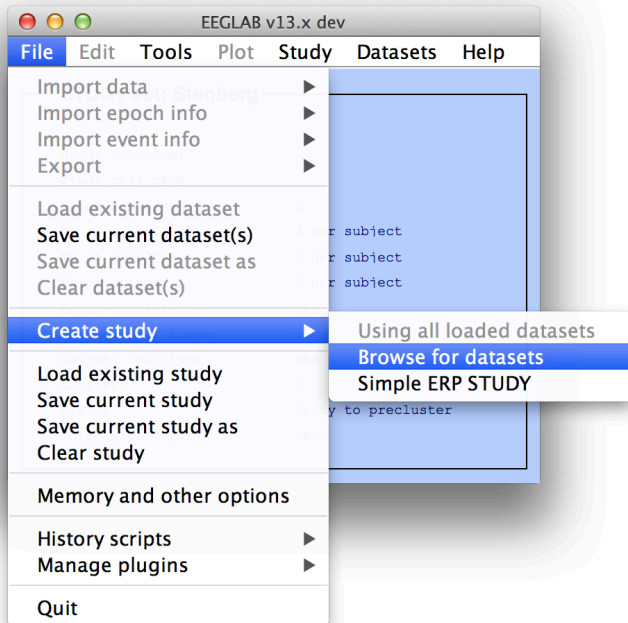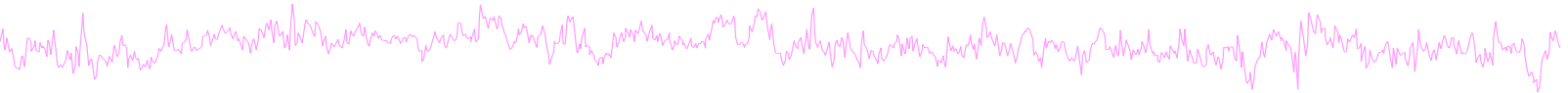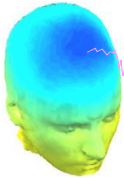

Create a new STUDY set -- pop_study()

Edit STUDY set information - remember to save changes

STUDY set name:

STUDY set task name:

STUDY set notes:

	dataset filename	browse	subject	session	condition	group	Select by r.v.	
1	/Users/arno/temp/STUDY/S01/	...	S01		Ignore		Comp.: 1 2 ...	Clear
2	/Users/arno/temp/STUDY/S01/	...	S01		Memorize		Comp.: 1 2 ...	Clear
3	/Users/arno/temp/STUDY/S01/	...	S01		Probe		Comp.: 1 2 ...	Clear
4		...						Clear
5		...						Clear
6		...						Clear
7		...						Clear
8		...						Clear
9		...						Clear
10		...						Clear

Important note: Removed datasets will not be saved before being deleted from EEGLAB memory

< Page 1 >


☐ Dataset info (condition, group, ...) differs from study info. [set] = Overwrite dataset info for each dataset on disk.

☐ Delete cluster information (to allow loading new datasets, set new components for clustering, etc.)

Help Cancel Ok


```
[STUDY ALLEEG] = std_editset( STUDY, ALLEEG, 'name','Stenberg','commands',{{'index' 1
'load' '/data/STUDY/S01/Ignore.set'} {'index' 2 'load' '/data/S01/Memorize.set'}
{'index' 3 'load' '/data/S01/Probe.set'} {'index' 1 'subject' 'S01'} {'index' 2
'subject' 'S01'} {'index' 3 'subject' 'S01'} {'index' 1 'condition' 'Ignore'} {'index' 2
'condition' 'Memorize'} {'index' 3 'condition' 'Probe'}},'updatedat','off' );
```

Looking at the function that create STUDY


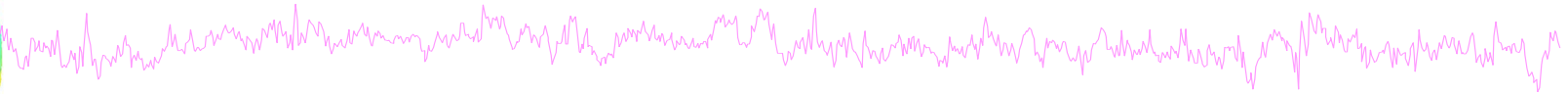


```
[STUDY ALLEEG] = std_editset( STUDY, ALLEEG, 'name','Stenberg','commands',  
    {'index' 1 'load' '/data/STUDY/S01/Ignore.set'} ...  
    {'index' 2 'load' '/data/S01/Memorize.set'} ...  
    {'index' 3 'load' '/data/S01/Probe.set'} ...  
    {'index' 1 'subject' 'S01'} ...  
    {'index' 2 'subject' 'S01'} ...  
    {'index' 3 'subject' 'S01'} ...  
    {'index' 1 'condition' 'Ignore'} ...  
    {'index' 2 'condition' 'Memorize'} ...  
    {'index' 3 'condition' 'Probe'}}, 'updatedat', 'off' );
```

```
[STUDY ALLEEG] = std_editset( STUDY, ALLEEG, 'name','Stenberg','commands',  
    {'index' 1 'load' '/data/STUDY/S01/Ignore.set' 'subject' 'S01' 'condition' 'Ignore'} ...  
    {'index' 2 'load' '/data/S01/Memorize.set' 'subject' 'S01' 'condition' 'Memorize'} ...  
    {'index' 3 'load' '/data/S01/Probe.set' 'subject' 'S01' 'condition' 'Probe'} ...  
    }, 'updatedat', 'off' );
```

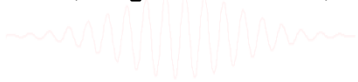

Exercise

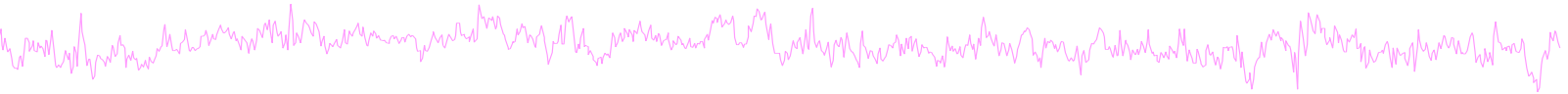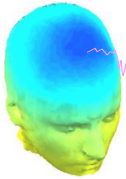

```
[STUDY ALLEEG] = std_editset( STUDY, ALLEEG, 'name', 'Stenberg', 'commands',  
    {{ 'index' 1 'load' '/data/STUDY/S01/Ignore.set' 'subject' 'S01' 'condition' 'Ignore' } ...  
    { 'index' 2 'load' '/data/S01/Memorize.set' 'subject' 'S01' 'condition' 'Memorize' } ...  
    { 'index' 3 'load' '/data/S01/Probe.set' 'subject' 'S01' 'condition' 'Probe' } ...  
    }, 'updatedat', 'off' );
```

1- Start EEGLAB and import the 3 datasets for Subject 1 (Ignore.set, Memorize.set and Probe.set) in a STUDY (menu Tools > Create STUDY > Browse for datasets)

2- Look in the history



3- Copy to a script, add also the “[ALLEEG EEG CURRENTSET ALLCOM] = eeglab;” at the beginning of your script and “eeglab redraw” after your command

4- Restart Matlab, execute the script, look at your STUDY info and design (menu *STUDY > Edit STUDY* info and *STUDY > Select/Edit STUDY design*)

5- Modify the script to import subject 1 to 4

6- Restart Matlab, execute the script, look at your STUDY info and design

Load dataset info from commandline

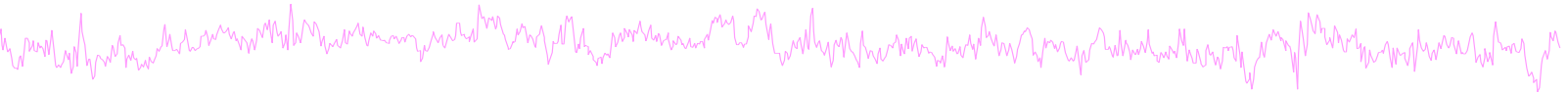
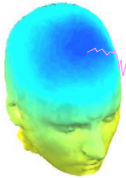


```
% Create Stern STUDY
[ALLEEG EEG CURRENTSET ALLCOM] = eeglab;
pop_editoptions( 'option_storedisk', 1);
subjects = {'S01' 'S02' 'S03' 'S04' 'S05' 'S06' 'S07' 'S08' 'S09' 'S10' 'S11' 'S12'};
filepath = '/Users/arno/temp/STUDY'; % XXXXX Change path here XXXXX
if ~exist(filepath), error('You need to change the path to the STUDY'); end;
commands = {}; % initialize STUDY dataset list


% Loop through all of the subjects in the study to create the dataset
for loopnum = 1:length(subjects) %for each subject
    IgnoreFile = fullfile(filepath, subjects{loopnum}, 'Ignore.set');
    MemorizeFile = fullfile(filepath, subjects{loopnum}, 'Memorize.set');
    ProbeFile = fullfile(filepath, subjects{loopnum}, 'Probe.set');
    commands = {commands{:} ...
        {'index' 3*loopnum-2 'load' IgnoreFile 'subject' subjects{loopnum} 'condition' 'Ignore'} ...
        {'index' 3*loopnum-1 'load' MemorizeFile 'subject' subjects{loopnum} 'condition' 'Memorize'} ...
        {'index' 3*loopnum 'load' ProbeFile 'subject' subjects{loopnum} 'condition' 'Probe'}};
end;
% Uncomment the line below to select ICA components with less than 15% residual variance
% commands = {commands{:} {'dipselect', 0.15}};
[STUDY, ALLEEG] = std_editset(STUDY, ALLEEG, 'name', 'Sternberg', 'commands', commands, 'updatedat', 'on');

% Update workspace variables and redraw EEGLAB
CURRENTSTUDY = 1; EEG = ALLEEG; CURRENTSET = [1:length(EEG)];
[STUDY, ALLEEG] = std_checkset(STUDY, ALLEEG);
eeglab redraw
```


STUDY structure

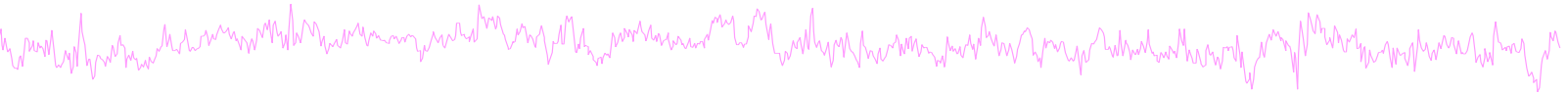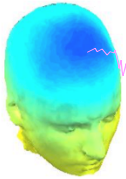

```
STUDY =  
      name: 'Sternberg'  
      task: 'Sternberg'  
datasetinfo: [1x39 struct]  
      notes: ''  
      filename: 'stern.study'  
      filepath: 'C:\Users\julie\Documents\Workshops\Finland\STUDY'  
      history: [1x7332 char]  
      subject: {1x13 cell}  
      group: {''}  
      session: []  
condition: {'ignore' 'memorize' 'probe'}  
      setind: [3x13 double]  
      etc: [1x1 struct]  
preclust: [1x1 struct]  
      cluster: [1x1 struct]  
changrp: [1x71 struct]  
      saved: 'yes'
```

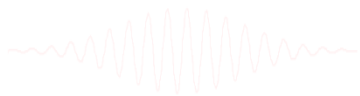

Understanding STUDY structure

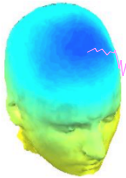

```
>> STUDY.datasetinfo(11) % access dataset 11  
ans =
```

```
    filepath: [1x61 char]
```

```
    filename: 'S04.set'
```

```
    subject: 'S04'
```

Subject 4!

```
    session: []
```

```
    condition: ''
```

```
    group: ''
```

```
    trialinfo: 1x350 struct
```

```
        index: 11
```

```
        comps: [1x24 double]
```

```
>> trialinfo(163) % access trial 163
```


```
ans =
```



```
    stimtype: 'Memorize'
```

```
    latency: 13201
```



```
    duration: 0
```

```
    ...
```

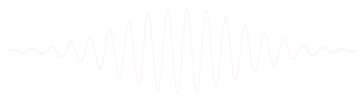

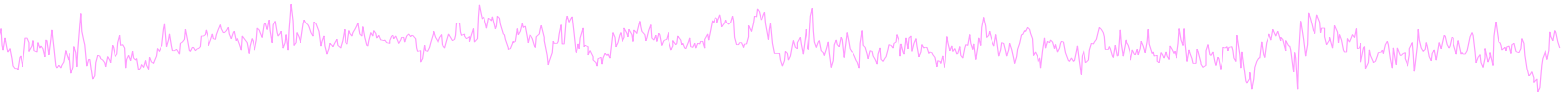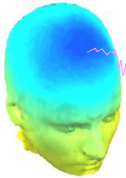


Select subjects

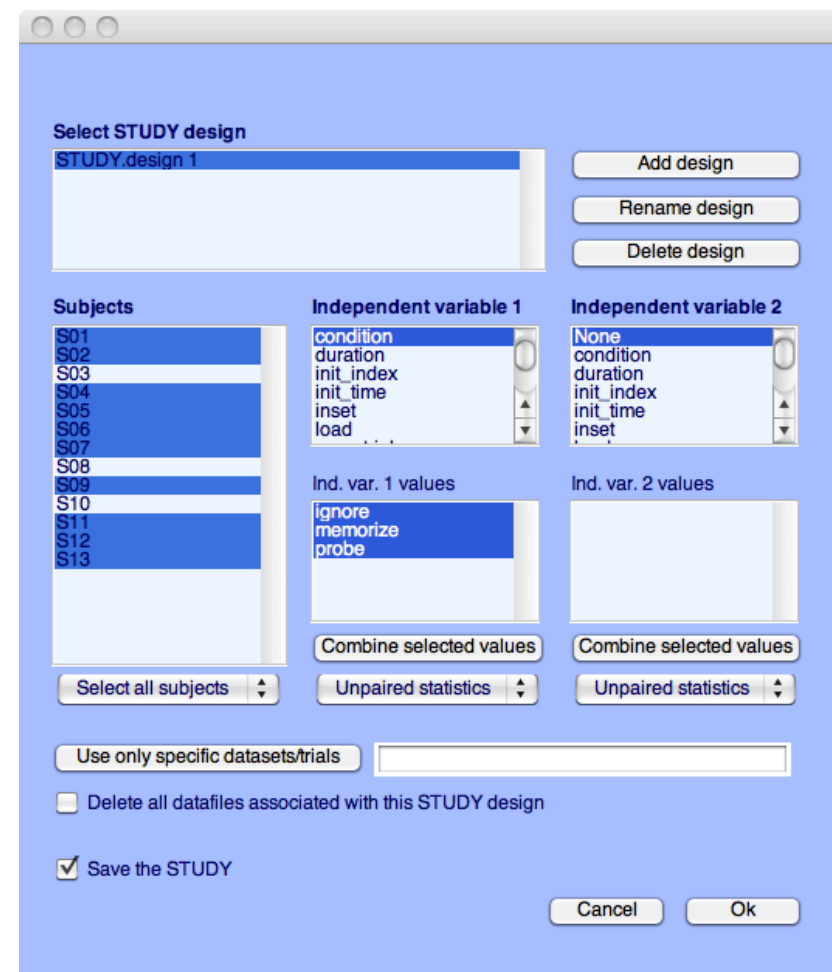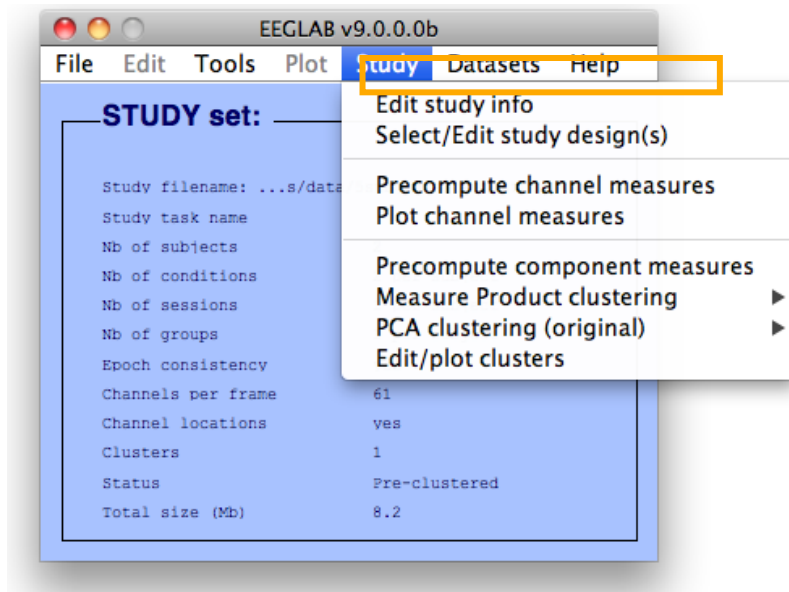

```
STUDY = std_makedesign(STUDY, ALLEEG, 3,
    'variable1','condition',
    'variable2','',
    'name','Design 3',
    'values1',{'ignore' 'memorize' 'probe'},
    'subjselect',{'S02' 'S03'},
    'dataselect',{'condition' {'probe'}});
```


STUDY design structure

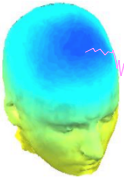

```
STUDY.design(1)

ans =

    name: 'Design 1 - compare letter types'
variable: [1x2 struct]
cases: [1x1 struct]
include: {}
cells: [1x39 struct]
```



Exploding the contents of each of these sub-structures, we obtain

```
    name: 'Design 1 - light and audio all subjects'

variable: [1x2 struct]
    (1).label : 'condition'
    (1).pairing: 'on'
    (1).value : {'ignore' 'memorize' 'probe'}
    (2).label : ''
    (2).pairing: 'off'
    (2).value : {}

cases: [1x1 struct]
    label: 'subject'
    value: {'S01' 'S02' 'S03' 'S04' 'S05' 'S06' }
```

Precompute data measures

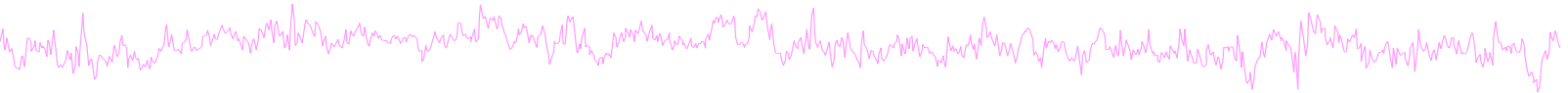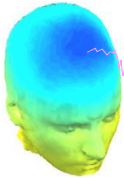

TIP: Compute all measures so you can test different combinations for clustering

Select and compute component measures for later clustering -- pop_precomp()

Pre-compute component measures for STUDY 'Sternberg'

☒ Compute ERP/spectrum/ERSP only for components selected by RV (set) or for all components (unset)


List of measures to precompute

☒ ERPs	Baseline ([min max] in ms)	[-200 0]	
☒ Power spectrum	Spectopo parameters		Test
☒ ERSPs	Time/freq. parameters	'cycles', [3 0.5], 'nfreqs', 100	Test
☒ ITCs			
☒ Scalp maps			

☐ Recompute even if present on disk

```
[STUDY ALLEEG] = std_precomp(STUDY, ALLEEG, 'components', 'erp', ...  
'on', 'rmbase', [-200 0], 'scalp', 'on', 'spec', 'on', ...  
'specparams', {}, 'ersp', 'on', 'erspparams', {'cycles' [3 0.5] ...  
'nfreqs', 100, 'freqs', [3 70]}, 'itc', 'on');
```

Precluster the data

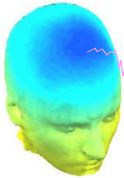

Select and compute component measures for later clustering -- pop_precl...

Build pre-clustering matrix for STUDY set: Sternberg

Select the cluster to refine by sub-clustering (any existing sub-hierarchy will be overwritten)

ParentCluster 1 (336 ICs)

Note: Only measures that have been precomputed may be used for clustering.

Measures	Dims.	Norm.	Rel. Wt.
☒ spectra	10	☒ 1	
☐ ERPs	10	☒ 1	
☒ dipoles	3	☒ 10	
☐ scalp maps	10	☒ 1	
☒ ERSPs	10	☒ 1	
☐ ITCs	10	☒ 1	
☐ Final dimensions	10		

Help

☐ Use Measure Product clustering





Help Cancel Ok

Additional settings for selected measures:

- Frequency range [Hz]: 3 25
- Time range [ms]:
- Use channel values: ☐ Absolute values
- Time range [ms]: 0 600
- Frequency range [Hz]: 3 30
- Time range [ms]:
- Frequency range [Hz]:



```
parentclust = 1;
[STUDY ALLEEG] = std_preclust(STUDY, ALLEEG, parentclust,{'spec','npca',5,...
    'norm',1,'weight',1,'freqrange',[3 25]},{ 'erp','npca',6,'norm',1,...
    'weight',1, 'timewindow',[0 400]},{ 'scalp','npca',10,'norm',1,'weight',1,...
    'abso',1},{ 'dipoles','norm',1,'weight',10},{ 'ersp','npca',20,...
    'freqrange',[3 30] , 'timewindow',[0 600], 'norm',1,'weight',1},{ 'itc',...
    'npca',6,'freqrange',[3 30], 'timewindow',[0 400] , 'norm',1, 'weight',1});
```

Cluster components

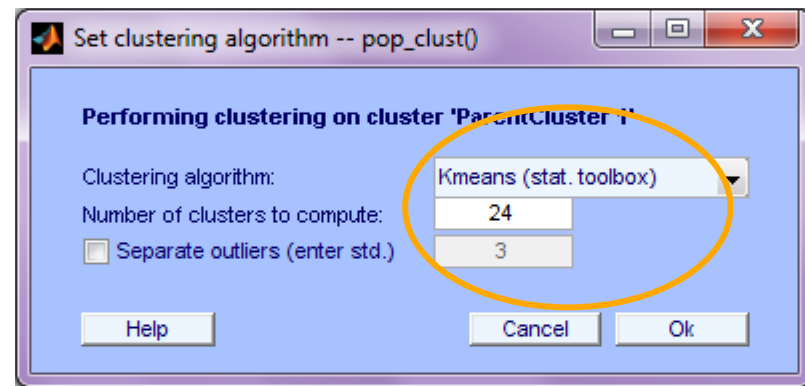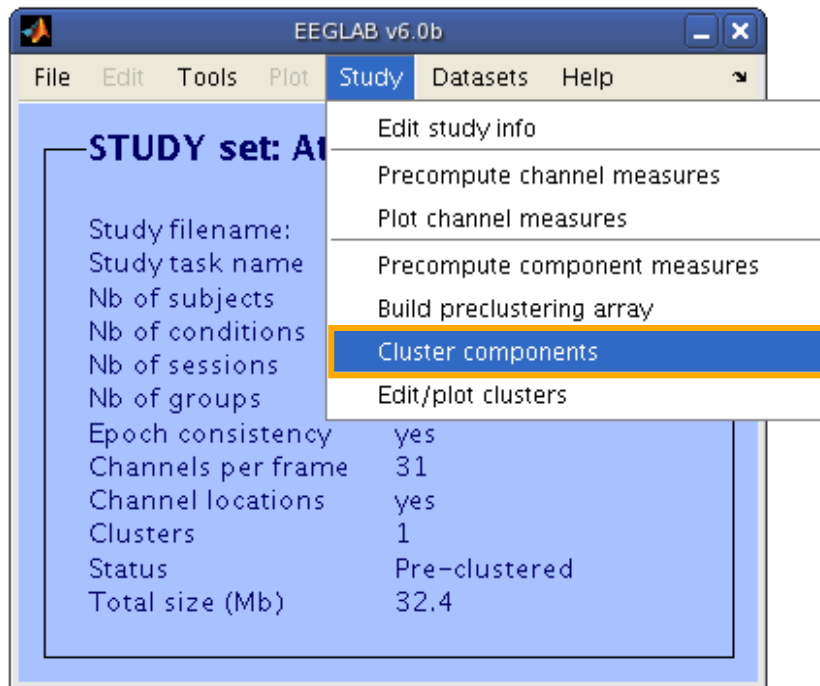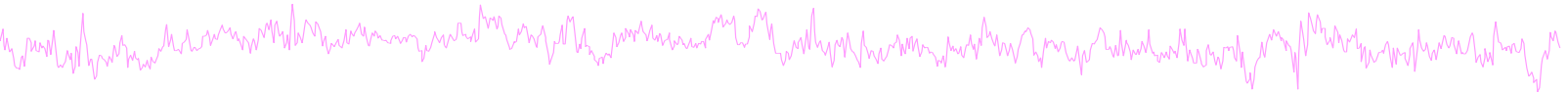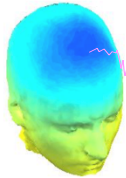

```
nclusts = 24; % choose # of clusters to create  
[STUDY] = pop_clust(STUDY, ALLEEG, 'algorithm', 'kmeans', 'clus_num', nclusts);
```

Understanding STUDY structure

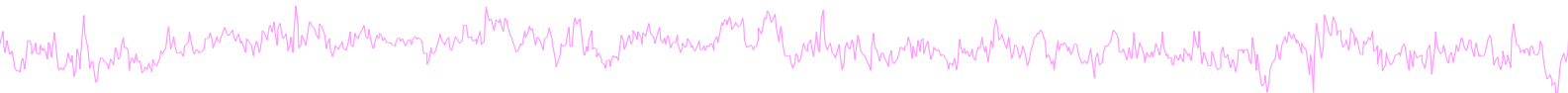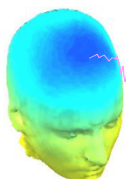

26 = # of clusters

```
>> STUDY.cluster  
1x26 struct array with fields:
```

```
parent  
name  
child  
comps  
sets  
algorithm  
preclust  
dipole  
allinds  
setinds
```

One cluster:

```
>> STUDY.cluster(6)  
ans =
```

```
parent: {'ParentCluster 1'}  
name: 'Cls 6'  
child: []  
comps: [35 7 12 35 10 23 7 30 4 ...]  
sets: [1 2 3 4 5 6 7 8 9 10 1 2 ...]  
algorithm: {'Kmeans' [24]}  
preclust: [1x1 struct]  
dipole: [1x1 struct]  
allinds: { [...] [...]; [...] [...]}  
setinds: { [...] [...]; [...] [...]}
```


IC indices

dataset indices for ICs

Design IC indices



Design indices

6 = cluster index

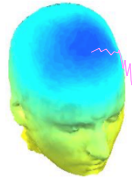


Choose which channel


Choose which subject

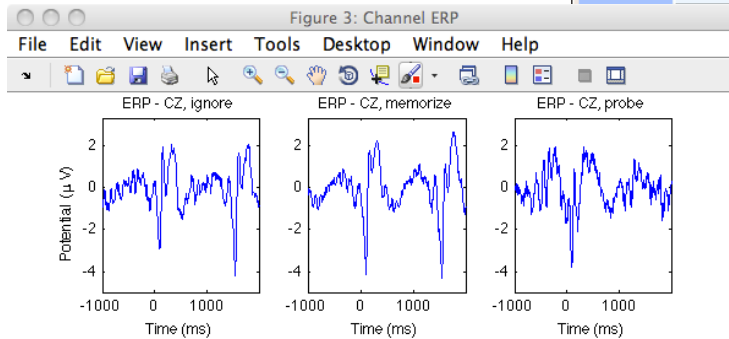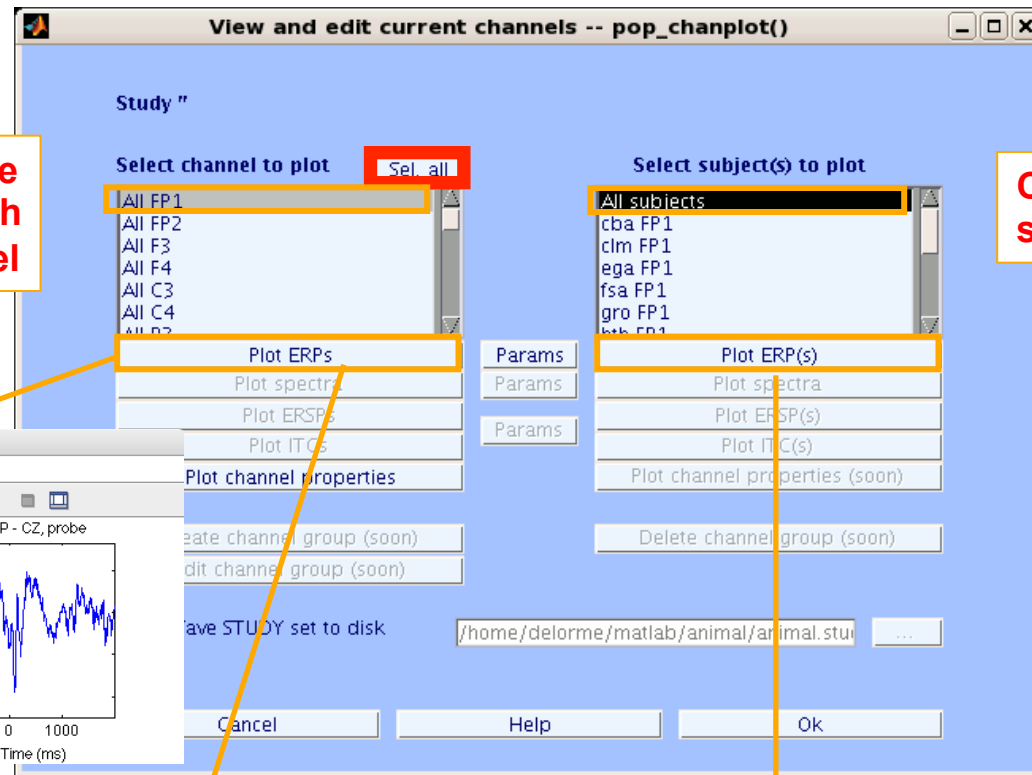

```
STUDY = std_erplot(STUDY,ALLEEG,'channels',{'FP1'});
```

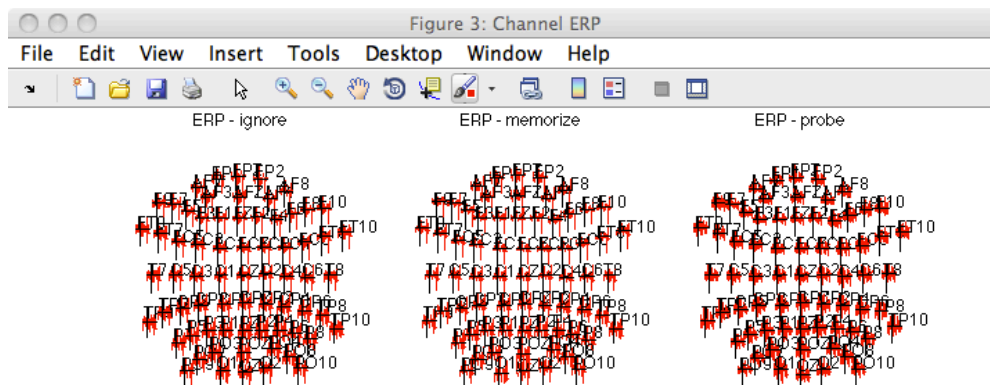


```
STUDY = std_erplot(STUDY,ALLEEG,'channels',{'FP1' ... });
```

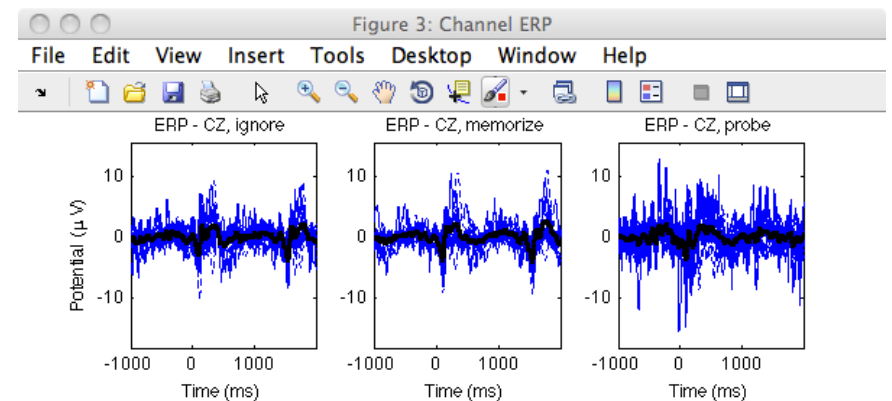



```
STUDY = std_erplot(STUDY,ALLEEG,'channels',{'FP1'}, 'plotsubjects', 'on' );
```

Advanced plotting features

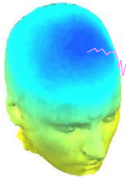

```
STUDY = std_erpplot(STUDY,ALLEEG,'channels',{ 'CZ' });  
std_plotcurve(STUDY.changrp(39).erptimes,  
STUDY.changrp(39).erpdata, 'plotconditions',  
'together', 'plotstderr', 'on', 'filter', 30);
```

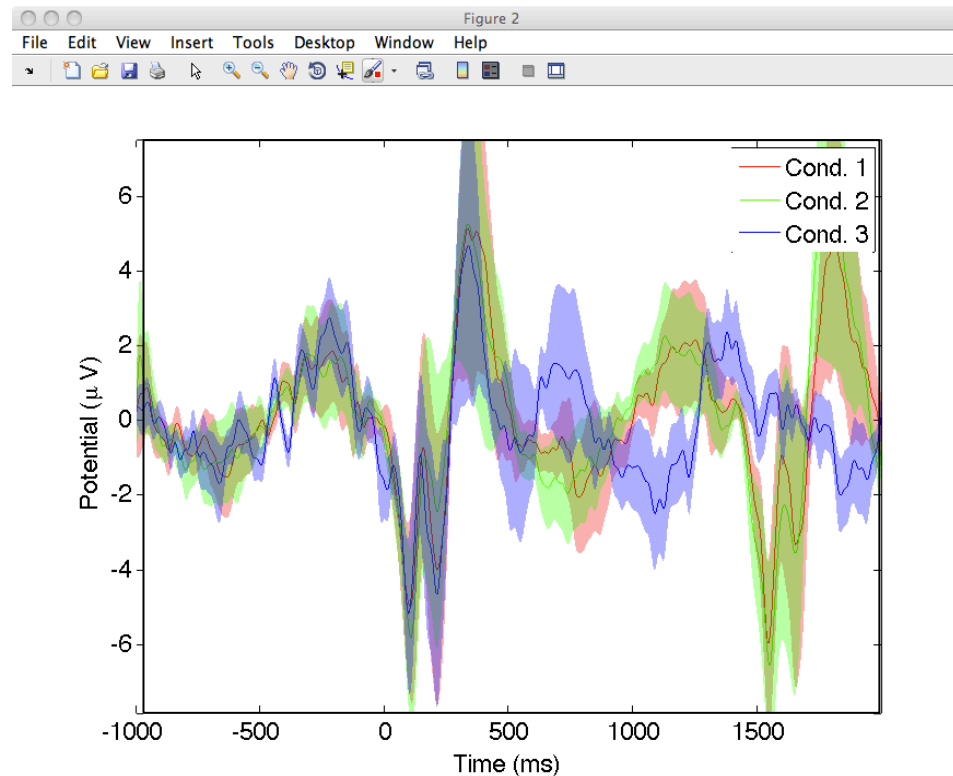

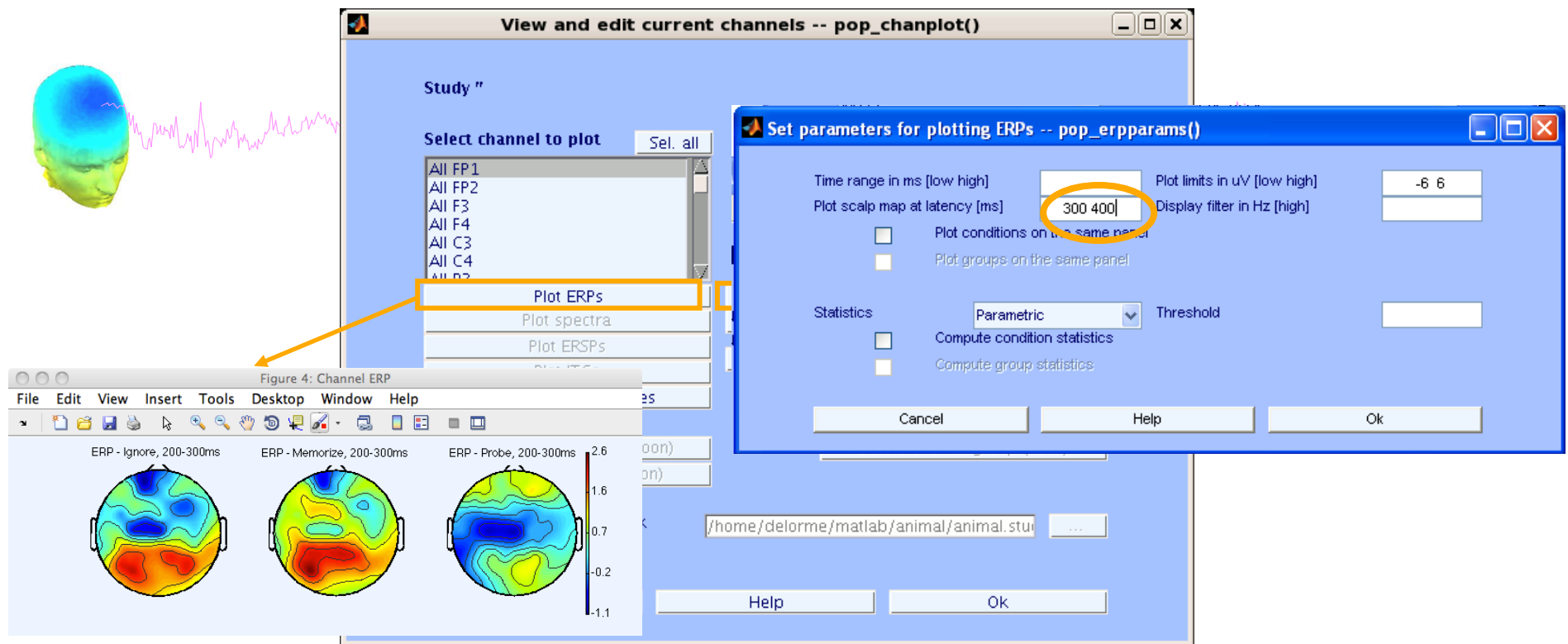

```
STUDY = std_erpplot(STUDY,ALLEEG, 'topotime',[200 300] , 'channels',{'OZ' 'O2' 'FP1' 'FPZ' 'FP2'});
[STUDY erpdata ] = std_erpplot(STUDY,ALLEEG, , 'topotime',[200 300] , 'channels',{'OZ' 'O2'});
```

Exporting to excell file

[1x67x13 double]

[1x67x13 double]


[1x67x13 double]

```
xlswrite('myxlsfile',squeeze(erpdata{1}),1);
```

```
xlswrite('myxlsfile',squeeze(erpdata{2}),2);
```

```
xlswrite('myxlsfile',squeeze(erpdata{3}),3);
```


Exporting text file

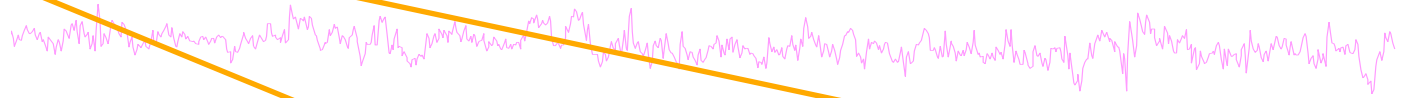

-0.13	-0.4	3.7	-0.9	-1.5	0.23	-0.98	1.8	2.3	-1.4	-2.8	-0.03	3.5
-0.54	-1.3	3.6	-1.1	-1.2	0.62	-0.91	1.6	2.2	-0.98	-7.7	-0.42	3.2
-0.77	-0.06	3.6	-1.4	-1.2	0.78	-0.91	1.2	2.1	-0.66	-0.76	-1	2.5
-0.61	-0.83	3.7	-1.2	-1.2	0.53	-0.88	1.1	1.7	-1.2	-1.8	-1.2	1.6
-0.34	-0.79	3.7	-0.98	-1.2	0.17	-0.72	1	1.4	-1.7	-2.3	-0.72	1.4
-0.27	-0.42	3.2	-0.69	-1.4	-0.04	-0.29	0.97	0.81	-2.5	-1.5	-0.38	1.7
0.097	-0.58	3.2	-0.61	-1.2	-0.32	0.36	0.47	2.1	-0.96	-2.8	0.89	2.4
0.43	-0.04	2.3	-0.47	-0.87	-0.37	0.21	0.83	3.1	-0.53	-0.85	1.2	3.4
0.21	-0.54	2.4	-0.07	-0.05	-0.08	-0.08	1	3.3	-0.42	-3.7	0.92	3.8
-0.1	-1.1	2.7	-0.33	-0.28	0.48	-0.5	1.2	3.3	-0.53	-2	0.36	4
-0.51	-2.2	2.9	-0.59	-0.23	1.3	-0.72	1.4	3.3	-0.14	-16	-0.05	3.9

```
dlmwrite('erpfile.txt',squeeze(erpdata{1}),'delimiter', '\t', 'precision', 2);
```

```
dlmwrite('erpfile.txt',squeeze(erpdata{2}),'-append', 'roffset', 1,  
'delimiter', '\t', 'precision', 2);
```

```
dlmwrite('erpfile.txt',squeeze(erpdata{2}),'-append', 'roffset', 1,  
'delimiter', '\t', 'precision', 2);
```

STUDY Script

```
% Create Stern STUDY
[ALLEEG EEG CURRENTSET ALLCOM] = eeglab;
pop_editoptions( 'option_storedisk', 1);
subjects = {'S01' 'S02' 'S03' 'S04' 'S05' 'S06' 'S07' 'S08' 'S09' 'S10' 'S11' 'S12'};
filepath = '/Users/arno/temp/STUDY'; % XXXXX Change path here XXXXX
if ~exist(filepath), error('You need to change the path to the STUDY'); end;
commands = {}; % initialize STUDY dataset list

% Loop through all of the subjects in the study to create the dataset
for loopnum = 1:length(subjects) %for each subject
    IgnoreFile = fullfile(filepath, subjects{loopnum}, 'Ignore.set');
    MemorizeFile = fullfile(filepath, subjects{loopnum}, 'Memorize.set');
    ProbeFile = fullfile(filepath, subjects{loopnum}, 'Probe.set');
    commands = {commands{:} ...
        ['index' 3*loopnum-2 'load' IgnoreFile 'subject' subjects{loopnum} 'condition' 'Ignore'] ...
        ['index' 3*loopnum-1 'load' MemorizeFile 'subject' subjects{loopnum} 'condition' 'Memorize'] ...
        ['index' 3*loopnum 'load' ProbeFile 'subject' subjects{loopnum} 'condition' 'Probe']};
end;



% Uncomment the line below to select ICA components with less than 15% residual variance
% commands = {commands{:} {'dipselect', 0.15}};
[STUDY, ALLEEG] = std_editset(STUDY, ALLEEG, 'name','Sternberg','commands',commands,'updatedat','on');

% Update workspace variables and redraw EEGLAB
CURRENTSTUDY = 1; EEG = ALLEEG; CURRENTSET = [1:length(EEG)];
[STUDY, ALLEEG] = std_checkset(STUDY, ALLEEG);
eeglab redraw

[STUDY ALLEEG] = std_precomp(STUDY, ALLEEG, {}, 'rmicacomps','on','interp','on','recompute','on','erp','on');
STUDY = pop_erpparams(STUDY, 'topotime',[200 300] );
[STUDY erpdata] = std_erpplot(STUDY,ALLEEG,'channels',{'LEYE' 'REYE' 'OZ' 'O2' 'FP1' 'FPZ' 'FP2' 'AF7' ...
    'AF3' 'AFZ' 'AF4' 'AF8' 'F9' 'F7' 'F5' 'F3' 'F1' 'FZ' 'F2' 'F4' 'F6' 'F8' 'F10' 'FT9' ...
    'FT7' 'FC5' 'FC3' 'FC1' 'FCZ' 'FC2' 'FC4' 'FC6' 'FT8' 'FT10' 'T7' 'C5' 'C3' 'C1' 'CZ' ...
    'C2' 'C4' 'C6' 'T8' 'TP9' 'TP7' 'CP5' 'CP3' 'CP1' 'CPZ' 'CP2' 'CP4' 'CP6' 'TP8' 'TP10' ...
    'P7' 'P5' 'P3' 'P1' 'PZ' 'P2' 'P4' 'P6' 'P8' 'PO9' 'PO7' 'PO3' 'POZ' 'PO4' 'PO8' 'PO10' 'O1'});

dlmwrite('erpfile.txt',squeeze(erpdata{1}), 'delimiter', '\t', 'precision', 2);
dlmwrite('erpfile.txt',squeeze(erpdata{2}), '-append', 'roffset', 1, 'delimiter', '\t', 'precision', 2);
dlmwrite('erpfile.txt',squeeze(erpdata{2}), '-append', 'roffset', 1, 'delimiter', '\t', 'precision', 2);
```


STUDY datafiles

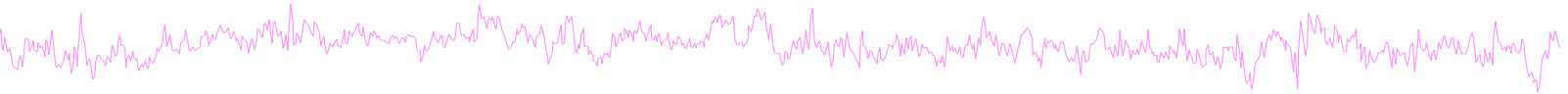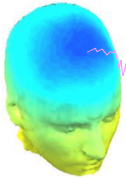



Datafiles

- *design1_S01_ignore.daterp* (ERPs)
- *design1_S01_ignore.datspec* (power spectra)
- *design1_S01_ignore.datersp* (ERSPs)
- *design1_S01_ignore.datitc* (ITCs)

- *design1_S01_ignore.icaerp* (ERPs)
- *design1_S01_ignore.icaspec* (power spectra)
- *design1_S01_ignore.icatopo* (scalp maps)
- *design1_S01_ignore.icaersp* (ERSPs)
- *design1_S01_ignore.icaitc* (ITCs)

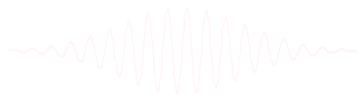


STUDY datafiles

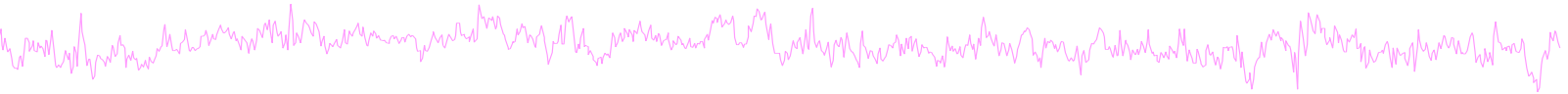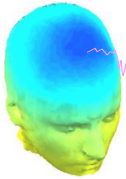



```
tmp = load('-mat', 'design1_S01_ignore.daterp');  
>> tmp
```

tmp =

```
      chan1:  [1x750 double]  
      chan2:  [1x750 double]  
      chan3:  [1x750 double]  
      chan4:  [1x750 double]  
      chan5:  [1x750 double]  
      chan6:  [1x750 double]  
      chan7:  [1x750 double]  
      ...  
      chan69: [1x750 double]  
      chan70: [1x750 double]  
      chan71: [1x750 double]  
      labels: {1x71 cell}  
      times:  [1x750 double]  
      datatype: 'ERP'  
      parameters: {'rmcomps' {1x1 cell} 'interp' [1x71 struct]}  
      datafile:  '/Volumes/donnees/data/STUDYstern/S01/Ignore.set'
```

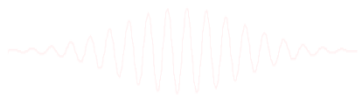


STUDY datafiles

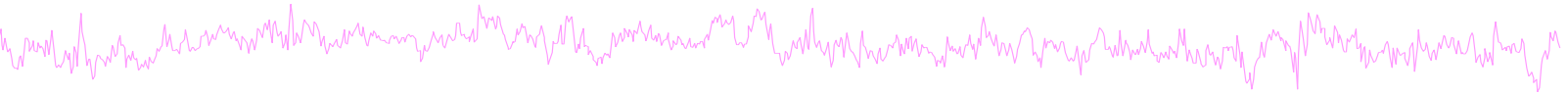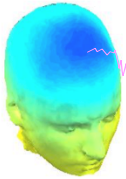



```
>> tmp = load('-mat', 'design1_S01_ignore.datersp');  
>> tmp
```

```
tmp =
```

```
      chan1_ersp: [50x100 single]  
chan1_erspbase: [1x50 single]  
chan1_erspboot: []  
      chan2_ersp: [50x100 single]  
chan2_erspbase: [1x50 single]  
chan2_erspboot: []  
      freqs: [1x50 double]  
      times: [1x100 double]  
      datatype: 'ERSP'  
parameters: {1x26 cell}  
chanlabels: {1x71 cell}  
datafile: '/Volumes/donnees/data/STUDYstern/S01/Ignore.set'
```

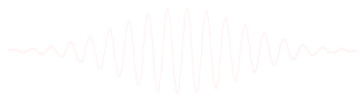


Exercises

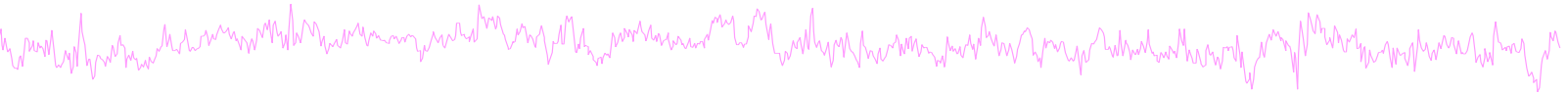
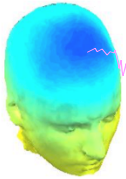


Suggestion for exercise


1. Load stern.study in STUDY folder
2. Choose a cluster to investigate and plot ERP
3. Get command line call for plotting ERP
4. Modify the script to plot spectrum instead of ERP
5. Export output to text file
6. Visualize output text file in another software

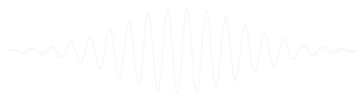



Advanced scripting: EEG pipeline

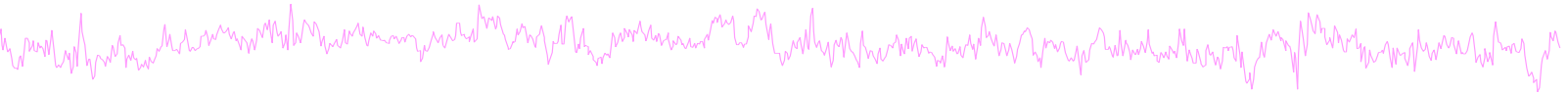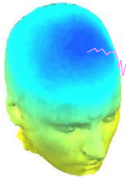

```
sInfo = [];  
sInfo(end+1).file = 'rawdata/S01.raw'; ← Raw data file  
sInfo(end).name = 'S01'; ← Subject name  
sInfo(end).bad_channels = { 'E1' }; ← Subject name
```

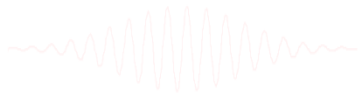


Advanced scripting: EEG pipeline

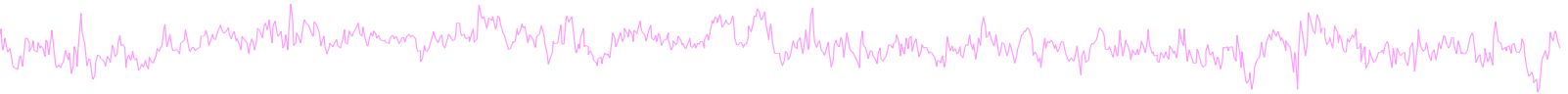
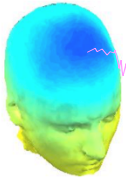


```
sInfo = [];  
sInfo(end+1).file = 'rawdata/S01.raw'; ← Raw data file  
sInfo(end).name = 'S01'; ← Subject name  
sInfo(end).bad_channels = { 'E1' }; ← Subject name  
sInfo(end).bad_data = [726 1495;6098 6831;13245 14057;15715 16399;22756 2445'
```

Copy the output from the eeg_eegrej function in the history

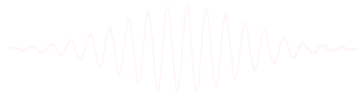


Advanced scripting: EEG pipeline

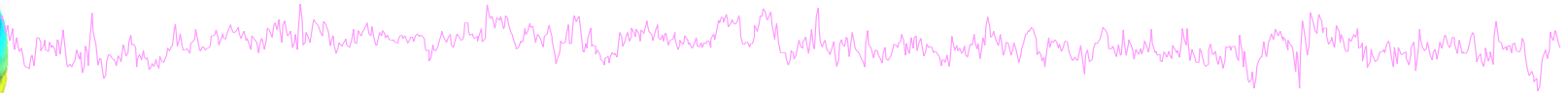
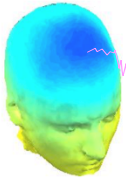


```
sInfo = [];  
sInfo(end+1).file = 'rawdata/S01.raw'; ← Raw data file  
sInfo(end).name = 'S01'; ← Subject name  
sInfo(end).bad_channels = { 'E1' }; ← Subject name  
sInfo(end).bad_data = [726 1495;6098 6831;13245 14057;15715 16399;22756 2445'  
sInfo(end).bad_comps = [1.6681 1.9870 0.3979 0.4444 -0.2274 -0.1433 -0.2626 .  
                        1.1917 -1.4838 0.7469 -1.1599 0.4773 -0.3257 0.3074 .
```

Copy transposed columns of the inverse weight matrix
EEG.icawinv for your selected artifact components

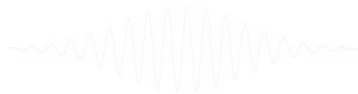


Advanced scripting: EEG pipeline

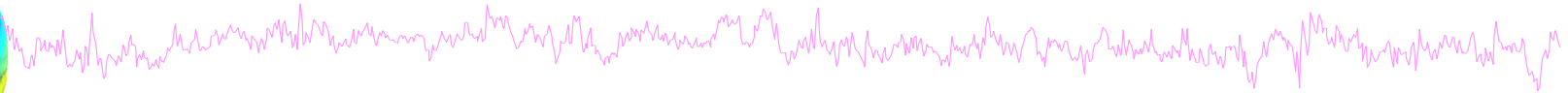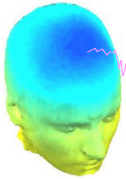

```
sInfo = [];  
sInfo(end+1).file = 'rawdata/S01.raw';  
sInfo(end).name = 'S01';  
sInfo(end).bad_channels = { 'E1' };  
sInfo(end).bad_data = [726 1495;6098 6831;13245 14057;15715 16399;22756 2445'  
sInfo(end).bad_comps = [1.6681 1.9870 0.3979 0.4444 -0.2274 -0.1433 -0.2626 .  
1.1917 -1.4838 0.7469 -1.1599 0.4773 -0.3257 0.3074 .  
  
sInfo(end+1).file = 'rawdata/S02.raw';  
sInfo(end).name = 'S02';  
sInfo(end).bad_channels = { };  
sInfo(end).bad_data = [41661 43713;24000 24833;44878 46501;48706 49210;51190  
sInfo(end).bad_comps = [0.6960 -0.8637 0.9087 -0.8028 0.4873 -0.2142 0.273'  
-0.0875 -0.4056 -0.0287 -0.3870 0.0600 -0.3716 0.342!  
2.1928 1.5712 0.8622 0.3215 -0.0357 -0.3125 -0.226!  
  
sInfo(end+1).file = 'rawdata/S03.raw';  
sInfo(end).name = 'S03';  
sInfo(end).bad_channels = { 'E10' 'E19' 'E20' 'E29' };  
sInfo(end).bad_data = [1 10449;19808 21815;25678 27254;29257 30010;34023 360:  
sInfo(end).bad_comps = [ 1.8583 2.0468 -0.0516 0.3159 -0.4256 -0.2770 -0.36'  
1.2189 -0.7385 1.2464 -0.8913 0.5475 -0.3971 0.29!  
-0.1248 -0.1358 -0.1954 -0.2533 -0.1555 -0.2313 -0.03!
```

datainfo.m file

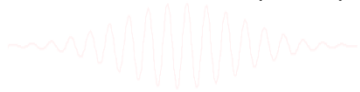

```

datainfo;
pop_editoptions( 'option_storedisk', 1);
Folder = 'ds_eeg_preprocessed_data';
if ~exist(outputfolder), mkdir(outputfolder); end;

for iSubject = 1:length(sInfo)

    % load dataset
    EEG = pop_biosig( sInfo(iSubject).file);
    EEG.setname = sInfo(iSubject).name;

    % preprocess data
    EEG = pop_iirfilt( EEG, 0.5, 0, [], 0, 0); % high pass filtering
    EEG = pop_iirfilt( EEG, 0, 55, [], 0, 0); % low pass filtering
    EEG = pop_select(EEG, 'nochannel', sInfo(iSubject).bad_channels); % remove bad channels
    EEG = pop_reref( EEG, []); % average reference (optional)
    EEG = eeg_eegrej( EEG, sInfo(iSubject).bad_data); % remove bad portions of data

    % run ICA
    EEG = pop_runica(EEG, 'pca', EEG.nbchan-1);


    % tag bad components
    if ~isempty(sInfo(iSubject).bad_comps)
        [corr,indx,indy] = matcorr(EEG.icawinv', sInfo(iSubject).bad_comps);
        if ~all(abs(corr(1:size(sInfo(iSubject).bad_comps,1))) > 0.92),
            error('Correlation too low'); end;
        EEG.reject.gcompreject(indx(1:size(sInfo(iSubject).bad_comps,1))) = 1;
    end;

    % extract data epochs
    EEG = pop_epoch(EEG, { 'S1' 'S2' });

    % save dataset
    EEG.saved = 'no';
    EEG = pop_saveset( EEG, 'filepath', folder, 'filename', [ sInfo(iSubject).name '.set' ]);
end;

```

Create STUDY

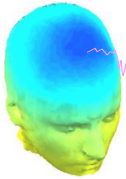

```
datainfo;
pop_editoptions( 'option_storedisk', 1);
outputEEGFolder = 'preprocessed_data';
studyCommand     = {};


% generate STUDY commands
for iSubject = 1:length(sInfo)
    fileName = fullfile(outputEEGFolder, [ sInfo(iSubject).name '.set' ]);
    studyCommand = [ studyCommand { 'index' iSubject 'load' fileName 'subject' ...
        sInfo(iSubject).name } ];
end;

% create data
[STUDY ALLEEG] = std_editset( [], [], 'name', 'test', 'commands', studyCommand, ...
    'updatedat','off', 'filename', 'test.study', 'resave', 'on');
STUDY = std_makedesign(STUDY, ALLEEG, 1, 'variable1','', 'variable2','', 'name', 'Design 1', ...
    'pairing1','on', 'pairing2','on', 'delfiles','off', 'defaultdesign','off', 'subjselect', ...
    {sInfo.name}, 'filepath', 'studyfiles');

% update workspace variables and redraw EEGLAB
CURRENTSTUDY = 1; EEG = ALLEEG; CURRENTSET = [1:length(EEG)];
[STUDY, ALLEEG] = std_checkset(STUDY, ALLEEG);
eeglab redraw

% precompute and plot data
[STUDY ALLEEG] = std_precomp(STUDY, ALLEEG, {}, 'interp', 'on', 'recompute', 'on', 'erp', 'on', ...
    'topotime', [200 300]);
[STUDY erp] = std_erpplot(STUDY, ALLEEG, 'channels', {allchanlocs.labels}, );
```

Exercise: build your own pipeline

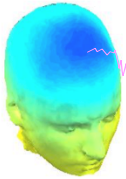

Suggestion for exercise

1. Load SimpleOddball.set dataset
2. High pass filter at 1Hz (menu Tools > Filter)
3. Reject bad portion of data by hand (menu Tools > Reject continuous...)
4. Re-reference to average ref. (optional) (menu Tools > Re-reference)
5. Build your *datainfo.m* file using EEGLAB history (see scripting lecture)
6. Run ICA (use PCA if necessary); select bad ICA components
7. Add components to you *datainfo.m* file
8. Create a STUDY with this single file
9. Compare the ERP for Oddball (type 2) and Standard (type 1) and use single-trial statistics with cluster correction for multiple comparisons
10. Build a script that creates the STUDY and perform the same analysis
11. Save the figure at the end of the script in eps or jpg format (“print –deps file” command or “print –djpg file” command).
12. Run the full pipeline (dataset processing and STUDY processing)
13. Change the filtering in the pipeline (step 2) and observe effects